


WIZLOC

Cervical Collar



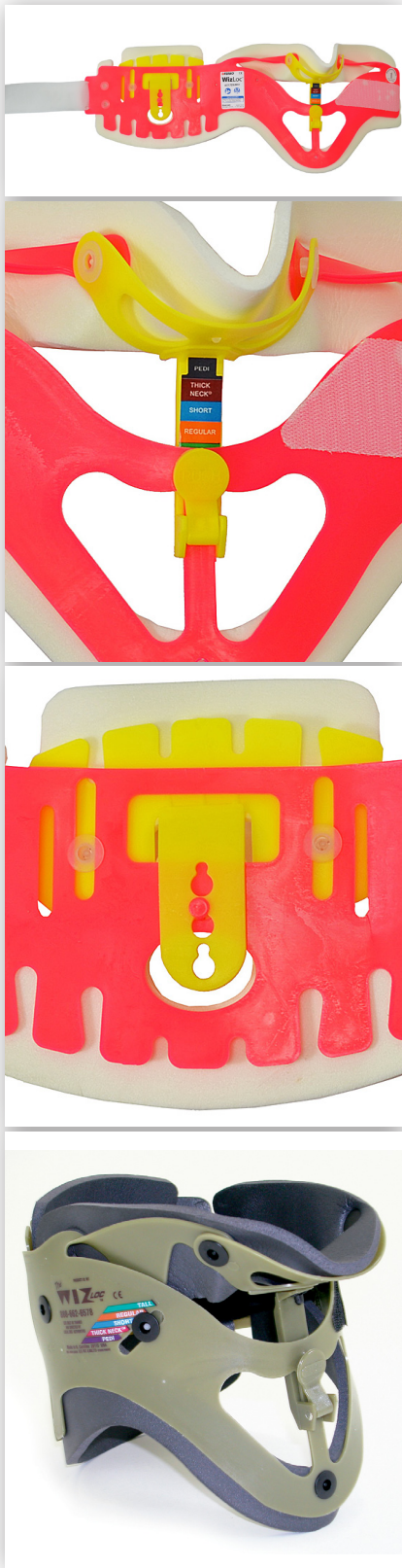

Fits 90% of all patients guaranteeing absolute stability of the cervical spine.

MAIN FEATURES

- The universal design of the WizLoc Cervical Collar fits a wide range of patient sizes ranging from 6 years old through adult.
- It features fast 3-step adjustment and application, with unique self-centering mandibular support and adjustable occipital support with keyhole lock pin. The occipital support is adjustable in 3 positions:
 - **Pedi-Tick Neck:** for pediatric patients, patients with a very short, thick neck, and patients with a low, well-defined occipital area;
 - **Regular-Short:** for most adult patients;
 - **Tall:** for patients with long necks and for patients with no defined occipital area.
- On-patient sizing for greater accuracy.
- High visibility colors for safety and identification.
- Equipped with Velcro fastenings.
- Provides minimal interference with X-ray or MRI procedures.
- Can be decontaminated for re-use up to 10 times.
- Available also in military version.

SPECIFICATION

Length max.	680 mm
Height (front)	170 mm
Height (rear)	200 mm
Weight	0,140 kg

MODELS

Description	Code
449-I cervical collar WizLoc (3 pcs)	0822074
449-I cervical collar WizLoc (military version)	0819759