

EASYFIX VACUUM MATTRESS

CODE: 271501001



EasyFix Vacuum Mattress



EasyFix Vacuum Mattress

EasyFix is a development work based on earlier experiences from the original high quality materials

EasyFix vacuum mattress make use of new body shape design to optimize the vacuum mattress capabilities.

EasyFix unique construction which deliver much improved fixation of upper back, neck and head by reinforced base plate that cover all of cervical spine. This together with fixation restraint system using x-system to support and create most effective immobilisation.

- Head, neck and upper back part improved with wood part to secure stronger immobilisation
- New improved chambers system
- Compliant to EN 1789 and 1865
- Carry capacity of 250 kg, 10 handles
- Light weighted and small storage

Specifications


Maximum shrinkage on length: 1%
 Maximum shrinkage on width: 3%
 Approved according to European Standard EN 1865 and EN 1789



X-ray and MRI transparent



Size:

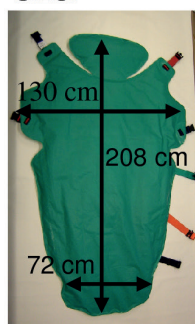

Weight:

8,0 kg without fixation straps
 8,5 kg with fixation straps
 9,6 kg total weight including fixation straps + footrest pump

Material

Outer material: Polyester fibre, PVC-coated (flame retardant)
 Inner bag: Polyester
 Granules: Styropor granules (flame retardant)
 Wood board upper part: Birch plywood

EasyFix	Dimensions and weight
EasyFix Vacuum Mattress (5 handles on each side)	L 208 x W130/72 x 4.5 cm weight 8,0 kg
EasyFix	Art. No.
Ordering specifications	
EasyFix Vacuum Mattress	271 501 000
EasyFix Vacuum Mattress Ready2go incl. Vacuum footrest pump+Repair Kit+ Bag	271 501 001

Features / Benefit	Easy Fix
Compliant to EN 1789 and 1865 Approved for European standard and crash tested for 10-G	X
Germa head pillow design including rigid wood part solution <i>Efficient, quick and improved neck and head immobilisation, with visible ears in one go</i>	X
One-piece back, <i>providing the patient with even spine support.</i> <i>Also allows preliminary X-ray diagnostics, direct in to x-ray no extra handling of patient</i>	X
Germa new chamber technology for quick and easy preparation <i>Gives rescuers quick handling and excellent immobilisation for the patient ensuring the granulate is evenly dispersed to give the patient stable and optimal immobilisation.</i>	X
Deliver consistent and rigid support and immobilisation of the patient <i>Safe, premium immobilisation and convenience of patient fixation.</i>	X
Handles ergonomic grip, <i>making patient lifts as comfortable as possible for rescuers.</i> Handles are reinforced so a <i>patient can be lifted in a patient lifting device.</i> Handles fixed to the bottom of mattress on piece, <i>maintain rigid immobilisation when lifting the patient.</i>	X
Robust welding—fully welded without any open seams or textiles, <i>making hygienic cleaning of the mattress easy.</i>	X
Valve with special designed hole, ensuring almost all suction devices fit directly into the valve, <i>no adapters are needed.</i> Reverse air lock system, <i>no closing of valve at stressful site.</i> Made from plastic makes it ideal for X-ray.	X
Improved V-shaped design <i>to make sure upper wide part will make all sized patients the very best immobilisation of the back, neck and head.</i>	X
X-designed restraint belt fixation close to the patient, <i>support immobilisation process and comfortable for patient.</i> Supplied with five colour coded straps either with fix lock, <i>helping the rescuer to place the straps correctly and easy to stretch and adjust.</i>	X
Outer material is made of flame retardant PVC-coated fabric. <i>It is flexible making it easy to mould the mattress around the patient.</i>	X
Design of the upper part allow the rescuer to place the arm(s) in two different options, <i>make sure to give the rescuer access to patient arm for treatment requested via patient arm(s)</i>	X
Extra layer built under the mattress, <i>safety, minimize risk of damage the fixation, last longer.</i> Robust and efficient material to hold vacuum for more than 24 hours, <i>safety for patient and make confiding immobilisation.</i> High load capacity 250 kg, <i>safe to handle heavy duty people</i>	X
Closing mechanisms of the straps, plastic buckle Fix lock, <i>allowing the user to quick tighten and adjust restraints on the patient.</i>	X
PVC coated fabric makes it easy to clean and disinfect with most common disinfectants, <i>easy to keep clean and ready for use.</i> Repair Kit easy to use for effective repair, <i>quick to repair no return needed</i>	X
Folds easily, 20% less storage space required compared to All In One Vacuum Mattress, <i>it is easy to fold and store the mattress in a protective bag and save storage space.</i>	X
High quality footrest aluminium vacuum pump evacuate with more than –700 hpa, <i>durable and quick and time saving. According to EN 1865</i>	X